

User Manual

This manual tells you how to use Full-SMS. It shows features, setup steps, instructions, and support data.

1. Introduction

Full-SMS helps you to analyze fluorescence measurements from HDF5 files. This manual tells you how to install, use, and maintain the software.

2. System Requirements

- **Web Browser:** Use Google Chrome or Mozilla Firefox.
- **Internet:** You must have a stable internet connection.
- **Data Format:** The software requires HDF5 measurement files.


3. Getting Access

1. Visit the Full-SMS website
2. Signup for a Full-SMS account
3. Login to your new account.


4. Get Started

This section tells you how to use Full-SMS after onboarding.


1. On the dashboard page locate the create workspace button.

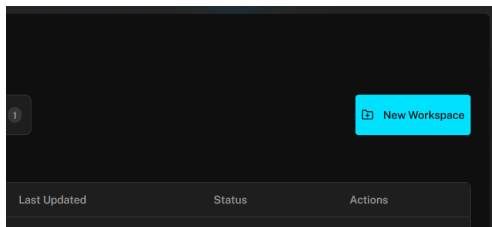

2. Follow the prompts on screen to create a new workspace.

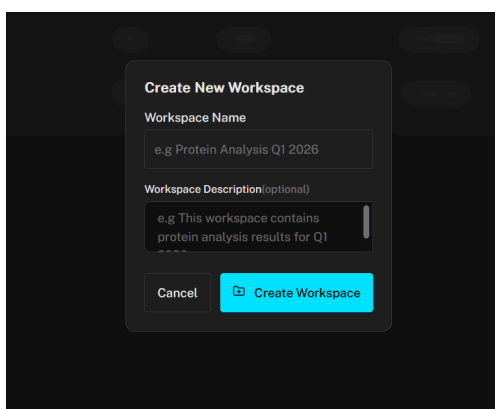

3. Your workspace will appear here shortly after creating it. Click on your workspace.

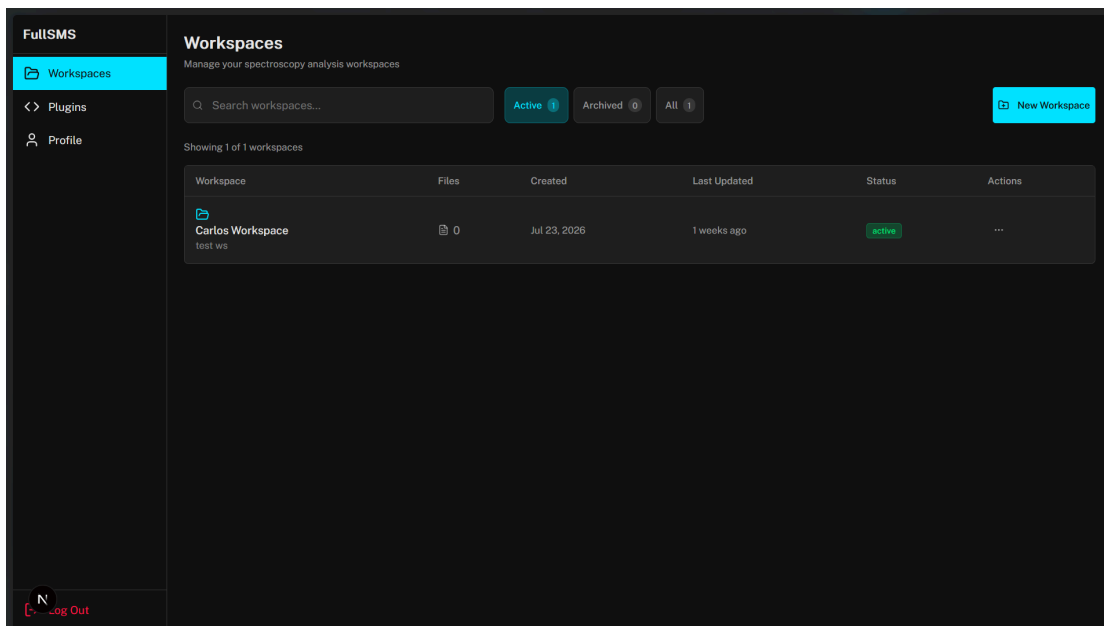


5. Use the Software


Read these steps to do the tasks:

Upload Data


1. Click the Upload button.

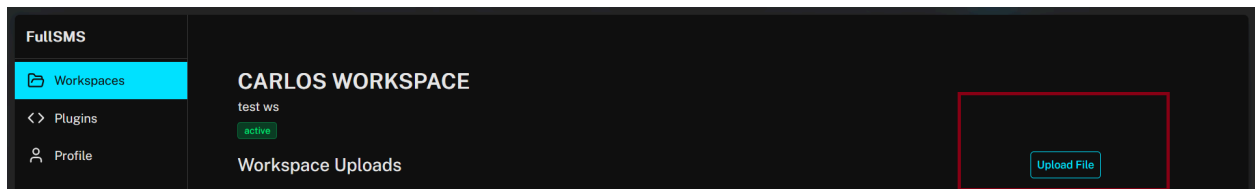

2. When the file upload modal opens, use it to select an .hdf5/.h5 file from your computer.

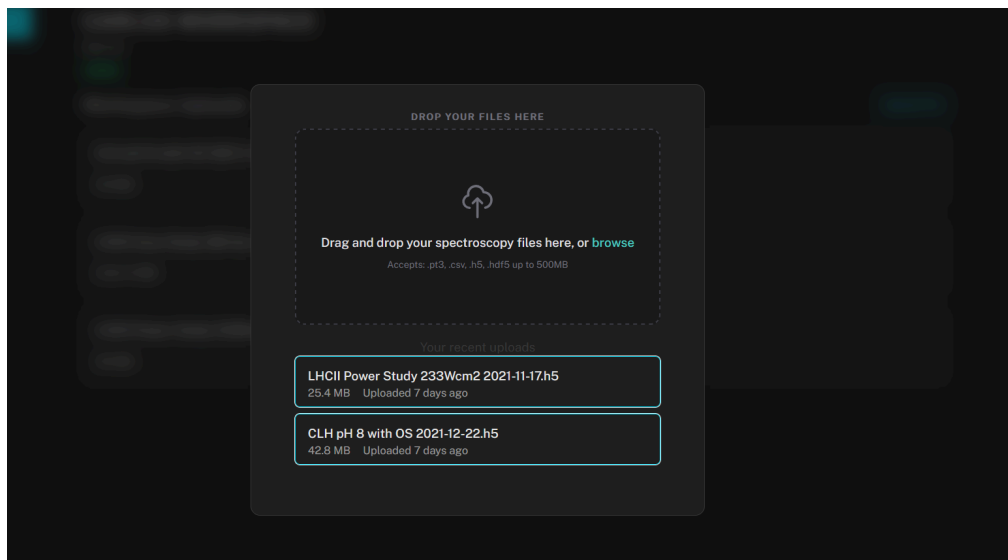


3. Click the Open button to get the data.

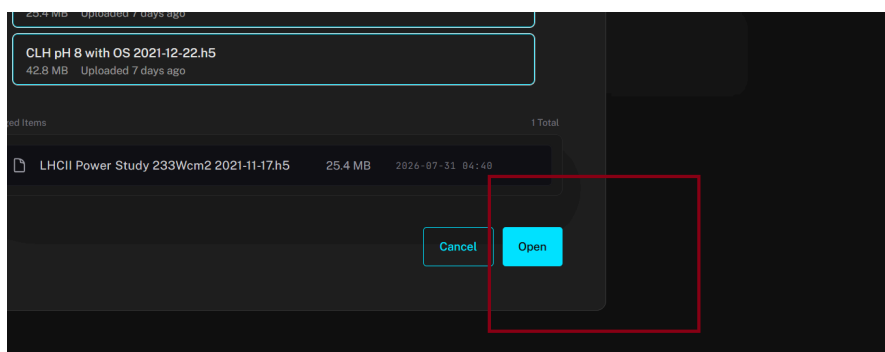


Run Change Point Analysis

1. When a file has been uploaded and parsed you will open it and navigate to the hub, there you will find the Intensity tab.

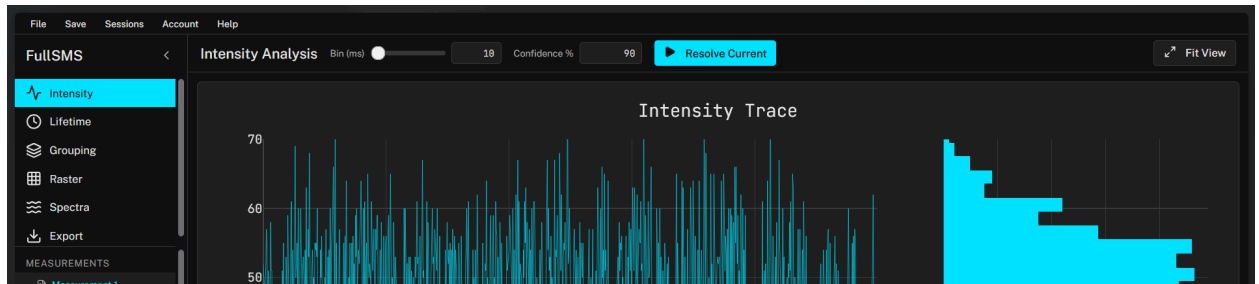

2. Click on the resolve current button.

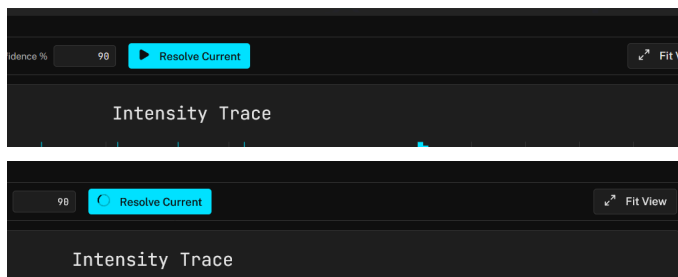

3. Once resolved, a red overlay will appear over your intensity trace.



6. Troubleshoot Issues

If you find a problem, read this table to find causes and solutions:

Problem	Cause	Solution
The file does not upload.	The file is not an HDF5 file.	Select an HDF5 file.
The lifetime fit fails.	The parameter values are incorrect.	Type the correct parameter values in the dialog.

7. Frequently Asked Questions (FAQ)

Can I use local scripts?

Technically, yes. Our plugin system allows you to inject your own custom analysis programs.

What data format does the software use?

The software uses HDF5/H5 format.

8. Document Changes

Version	Date	Changes
1.0	2026-07-31	Initial document release.